




SHIELDskin XTREME™
A REVOLUTION IN GLOVE TECHNOLOGY

DI++

EXTREME
CONTAMINATION CONTROL

SHIELDskin XTREME™

White Nitrile 300 DI++

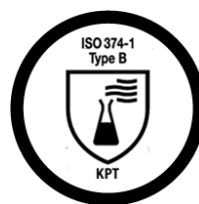



DI++

Extreme
Contamination
Control

- ⇒ Powder-free multiple DI washed ambidextrous standard length (300 mm / 11.8") non-sterile nitrile cleanroom gloves.
- ⇒ Personal Protective Equipment Category III (PPE - Complex Design) according to Regulation (EU) 2016/425.
- ⇒ Fully compliant to the latest PPE Protective gloves EU norms against chemicals, micro-organisms and viruses.

DESCRIPTION							
Formulation	Nitrile synthetic rubber (<i>acrylonitrile butadiene</i>).						
Design	White, ambidextrous, beaded cuff, textured fingertips.						
Packaging	100 gloves per PE bag - 10 bags per polybag - 1 polybag per carton.						
SIZES		6/XS	7/S	8/M	9/L	10/XL	11/XXL
Codes		69 8851	69 8852	69 8853	69 8854	69 8855	69 8856

STANDARDS	
CE registration	PPE Category III (Complex Design) - Regulation (EU) 2016/425. Notified Body No 0598: SGS Fimko Oy, Helsinki - FINLAND.
EU PPE norms	ISO 21420:2020, ISO 374-1:2016+A1:2018, ISO 374-2:2019, ISO 374-4:2019, ISO 374-5:2016, EN 16523-1:2015+A1:2018 and ISO 16604:2004 Procedure B.
EU MDR norms ¹	EN 455-1:2020, EN 455-2:2015, EN 455-3:2015 and EN 455-4:2009.
USA standards	ASTM D3767-03 (2020), ASTM D573-04 (2019), ASTM D412-16 (2016) and IEST-RP-CC005.4 (2013).
Other standards	EN 1149-1/2/3 & 5, ISO 10993-10:2010.

¹With reference to Regulation (EU) 2017/425 for Medical Devices

QUALITY	
Quality assurance	Production management in accordance with ISO 9001:2015 and ISO 13485:2016. Environmental management systems in accordance with ISO 14001:2015.
Technology	uniSHIELD™ single-walled protection to offer an ideal compromise between comfort and protection. Compatible with clean processing environments due to paperless packaging and multiple post leaching of gloves (multi washed in deionised water).

DOCUMENTATION	
Declaration of conformity	These documents can be freely downloaded from the product page on our website: www.shieldscientific.com .
EU type examination certificate	For easy access, scan the QR code.
User's instructions	
Certificate of conformance	To access CoC, you need to be registered. Please contact us at info@shieldscientific.com or call your SHIELD Scientific representative.



PHYSICAL PROPERTIES



NOMINAL THICKNESS	mm ²	mil	Norm
⇒ Finger	0.19	7.5	ASTM D3767-03 (2020)
⇒ Palm	0.14	5.5	
⇒ Cuff	0.11	4.3	

² Thickness (+/- 0.03 mm)

LENGTH	Minimum	Typical	Norm
⇒ From middle finger tip to edge of cuff	≥ 295 mm / 11.6"	305 mm / 12.0"	ISO 21420:2020

STRENGTH PROPERTIES	Force at break (spec.)		Ultimate elongation (spec.)	Force at break (typical)	Norm
⇒ Before aging	≥ 6.0N	14 MPa	≥ 500%	11.0N	EN 455-2:2015 ASTM D573-04 (2019) & ASTM D412-16
⇒ After aging	≥ 6.0N	14 MPa	≥ 400%	10.0N	

FREEDOM FROM HOLES	Performance	Norm
⇒ Acceptable Quality Level (AQL)	< 1.5 ³ - Level 2	ISO 374-2:2019

³ AQL as defined per ISO 2859-1:1999 for sampling by attributes.

RISKS	Description	Norm
Micro-organisms	1000 ml water test. Performance level 2, AQL < 1.5 (inspection level G1).	ISO 374-2:2019
Viruses	Viral penetration test using Phi-X174 bacteriophage according to ISO 16604:2004 Procedure B.	ISO 374-5:2016
Chemicals	<u>Performance</u> : Type B (KPT). <u>Permeation</u> : Extensively tested. Online chemical resistance guide on www.shieldscientific.com . <u>Degradation</u> : Tested for determination of resistance to degradation by chemicals.	ISO 374-1:2016+A1:2018 EN 16523-1:2015+A1:2018 ISO 374-4:2019

CLEANLINESS PROPERTIES

PARTICLES	Specification	Typical value	Test method
Particles/cm ² ≥ 0.5µm	< 850 particles	600 particles	IEST-RP-CC005.4

EXTRACTABLES (ION)	Specification (µg/cm ²)	Typical value (µg/cm ²)	Test method
Ammonium (NH ₄)	0.020	0.009	IEST-RP-CC005.4
Bromide (Br)	0.020	< 0.008	
Calcium (Ca)	0.200	0.120	
Chloride (Cl)	0.100	0.070	
Fluoride (F)	0.010	< 0.008	
Magnesium (Mg)	0.010	< 0.008	
Nitrate (NO ₃)	0.100	0.070	
Nitrite (NO ₂)	0.020	< 0.008	
Phosphate (PO ₄)	0.020	< 0.008	
Potassium (K)	0.030	0.013	
Sodium (Na)	0.020	0.011	
Sulphate (SO ₄)	0.020	< 0.008	
Zinc (Zn)	0.050	0.030	

EXTRA TESTS	Description	Test method
NVR	Maximum 30 µg/g.	IEST-RP-CC005.4
FTIR	Silicone free and non-detectable levels of amide and DOP.	IEST-RP-CC005.4
ESD	Tested for electrostatic properties.	EN 1149-1/2/3 & 5

ALLERGIES	
Bio-Compatibility	Demonstrated by skin irritation and sensitization tests in accordance with ISO 10993-10:2010.
Accelerators	Accelerator-free to minimize the risk of allergic contact dermatitis (also known as Type IV, delayed hypersensitivity or chemical allergy).
Chemical Allergens	Non-detectable levels using aqueous solution extraction (Phosphate buffered solution) and High Performance Liquid Chromatography (HPLC) assay method for quantitative analysis.
Latex Protein	Latex-free.